CONTADOR DIGITAL LEUCOCITARIO DIGITAL LEUKOCYTE COUNTER

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This manual should be available for all users of these equipments. To get the best results and a higher duration of this equipment it is advisable to read carefully this manual and follow the processes of use.

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AUXILAB S.L.

DECLARATION OF CONFORMITY
DIGITAL LEUKOCYTE COUNTER of Auxilab, S.L for the Directive
of Machines (89/392/CEE modified) and the regulations
adopted for their transposition.

NAME OF THE MANUFACTURER / IMPORTER:

AUXILAB, S.L.

ADDRESS:

Polígono Morea Norte, 8
31191 Beriain (Navarra)

WE STATE THAT:

DIGITAL LEUKOCYTE COUNTER
CODE 36048100

Is designed and manufactured according to:

- ♦ Directive 89/392/CEE, including the modifications and national regulations that transpose them.
- ♦ Directive 89/336/CEE modified over the electromagnetic compatibility.
- ♦ Directive 73/23/CEE modified over the electric security.

And that the following harmonized rules have been applied (or a part of them):

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BERIAIN 14th March 2007

Signed by: ALFONSO AINCIBURU SANZ
 DIRECTOR/MANAGER

Polígono Morea Norte, 8 31191 Beriain (Navarra) - Spain. Tel. 948 310 513 Fax 948 312 071
Internet: www.auxilab.es · Email: correo@auxilab.es

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- ◆ Always use original components and supplies. Other devices can be similar but they can damage the equipment.
- ◆ The leukocyte counter is supplied with a wire. It has to be plugged to an earth connection and the socket should be handy and ready to unplug the equipment in case of emergency.
- ◆ Never try to repair the leukocyte counter by yourself, since you will lose the warranty and may provoke damages to the general operating system or the electrical installation, as well as injuries to the people that usually handle the equipment (burns, hurts...) or damages in nearby equipments.
- ◆ In the event of breakdown please contact your distributor to overhaul through Nahita Technical Assistance Department.

Cleaning

- ◆ Never use scourers or substances that can grate for cleaning metallic parts such as stainless steel, aluminium, coatings, etc. as they damage the leukocyte counter and produce an early ageing of the equipment.
- ◆ Use a fluff-free cloth dampened with soaped water that does not contain abrasives.

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% key is pressed again, percentage decimals will be displayed.

For example, if a cell type represents 12.63 % of the total count, by pressing once the % key the number 12 will be displayed and by pressing it again 63 will be displayed.

4.To perform a new count, press RESET.

5.When operation is finished, turn the equipment off by pressing the ON/OFF key and disconnect it from the net.

NOTE: If you are not using the leukocyte counter for a long period of time please make sure it is disconnected from the net and protected from dust (this way you will avoid accidents and will extend its working-life).

Security

- ◆ The leukocyte counter must be used by previously qualified staff that knows how the equipment works thanks to the user manual.
- ◆ Put the leukocyte counter on top of a horizontal, plane and stable table making a free space at least at 30 cm per side.
- ◆ Do not put the leukocyte counter near any warm supply (burners, blowlamps...), nor expose it directly to the sun, etc. Avoid vibrations, dust and dry environments.
- ◆ When you are not using the leukocyte counter for a long period of time please make sure it is unplugged in order to avoid possible accidents.
- ◆ During operation, dangerous material such as flammable or infectious liquids must be out of the working area.
- ◆ It is essential to have the equipment switched off and unplugged from the net before cleaning, checking components or replacing any piece (e.g. replacement of a fuse).
- ◆ Never try to repair the leukocyte counter by yourself, since you will lose the warranty and may provoke damages to the general operating system or the electrical installation, as well as injuries to the people that usually handle the equipment (burns, hurts...).
- ◆ Try not to spill any liquid on the control panel, though it is properly insulated. In case you have any doubts do please immediately unplug the equipment (see Maintenance).
- ◆ Made under the European regulations for electrical security, electromagnetic compatibility and security on machines.

5. MAINTENACE AND CLEANING

To get the best results and a higher duration of the leukocyte counter it is essential to follow the processes of use.

Note: All the processes of use mentioned below will not have any value unless you keep a continued and careful maintenance.

- ◆ Please follow the processes of use of this manual.
- ◆ This manual should be available for all users of this equipment.
- ◆ No use objetos punzantes como bolígrafos etc, para tocar las teclas del panel, use sola

- ◆ We will only accept any equipment return within 15 days after delivery and provided it comes in its original wrapping.

Installation

- ◆ Before using the pH-meter, it is convenient for you to familiarize with its components and basic essentials, as well as with its control functions.

PLEASE READ THOROUGHLY THE INSTRUCTIONS BEFORE CONNECTING AND OPERATING WITH THIS EQUIPMENT WITH THE AIM OF ACHIEVING THE HIGHEST FEATURES AND THE MAXIMUM DURABILITY OF THE INSTRUMENT.

- ◆ Please put the leucocyte counter on top of a horizontal, plane and stable table making a free space at least at 30 cm per side. Do not put the leukocyte counter near any warm supply (burners, blowlamps...), nor expose it directly to the sun, etc.
- ◆ During operation, dangerous material such as flammable or infectious liquids must be out of the working area.
- ◆ Please insert the wire that feeds the AC electric current in the base of current 220 V, 50 Hz $\pm$ 10% provided with earth wire.
- ◆ Neither the manufacturer nor the distributor will assume any responsibility for the damages produced to the equipment during its installation or damages to people suffered by the improper use of the electric connection. The tension should be 220 V, 50 Hz $\pm$ 10%.

1. Turn the equipment on by pressing the ON/OFF switch. All displays will show "0".

NOTE: If 0 is not shown in all displays, press RESET.

2. Press the corresponding key for each leukocyte type (see chart).

With each keystroke a beep sounds and a unit is added to its corresponding display and also to the total count display.

NOTE: A maximum of 100 cells can be counted for each cell type, with a total maximum of 1000 cells.

Key	Leukocyte type	Key	Leukocyte type
1	Basophils	7	Promyelocytes
2	Variety of lymphocytes	8	Segmented neutrophils (mature)
3	Eosinophils	9	Ortocratic normoblasts
4	Monocytes	--	Other
5	Stab neutrophils (immature)		
6	Lymphocytes		

When the maximum total count (1000) is reached a continuous beep will sound.

3. Once count is finished, percentage of each cell type regarding the total count can be calculated. To do so, press % key and the percentage units and tens will be displayed.

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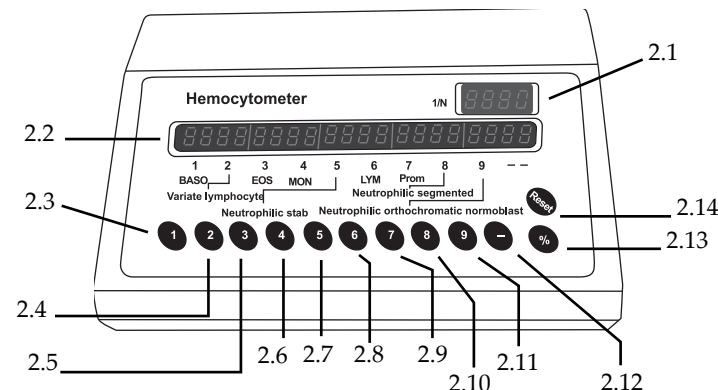
3. TECHNICAL SPECIFICATIONS

Leukocyte types	10
Count	0-100 for each cell type
Total count	1000
% count	Yes

4. INSTALLATION / SETTING UP

Preliminary inspection

Unwrap the leukocyte counter, take off the involving plastic and take off the polispan protection in which it comes fitted. Take off all protections and without connecting the leukocyte counter to the power supply, make sure that the equipment does not present any damage because of the shipment. In case it presents any damage, tell it immediately to your transport agent or dealer so that they can make the claims in the correct time limit.

Please keep the original wrapping; you will always need it for returns enclosed with all the accessories supplied.

Please check that all the accessories are enclosed with the equipment:

-Digital leukocyte counter unit.

INDEX OF LANGUAGES

Spanish	2-8
English.....	9-15

INDEX OF CONTENTS


1. USES OF THE INSTRUMENT	10
2. DESCRIPTION.....	10
3. TECHNICAL SPECIFICATIONS	11
4. INSTALLATION / SETTING UP.....	11
5. MAINTENANCE AND CLEANING.....	13
ANNEX I: CE CERTIFICATE.....	15

1. USES OF THE INSTRUMENT

Leukocyte count in blood samples is essential to determine anomalies that indicate pathological processes such as infections or leukaemia. To facilitate counting we offer a useful digital counter equipped with 10 operating keys for different types of leukocytes (basophils, eosinophils, monocytes, lymphocytes, mature and immature neutrophils, etc.). Pressing each key adds one unit to its corresponding display and to the total record. The counter also calculates the percentage for each cell type respecting to the total count and is equipped with a reset function.

The leukocyte counter is suitable for using in research, university and hospital laboratories.

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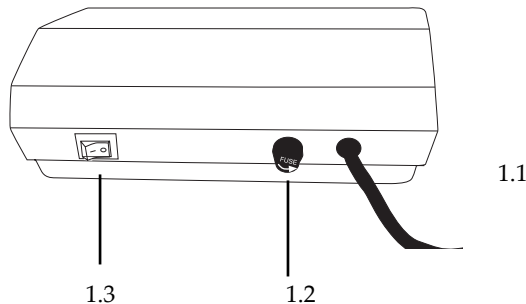Thank you for choosing this equipment. We sincerely wish that you enjoy your Nahita digital leukocyte counter. We highly recommend looking after this equipment according to what is stated in this manual.

Nahita develops its products according to the CE marking regulations as well as emphasizing the ergonomics and security for its user.

The correct using of the equipment and its good quality will permit you to enjoy this equipment for years.

The improper use of the equipment can cause accidents and electric discharges, circuit breakers, fires, damages, etc. Please read the point of Maintenance, where we expose the security notes.

TO GET THE BEST RESULTS AND A HIGHER DURATION OF THE EQUIPMENT IT IS ADVISABLE TO READ THOROUGHLY THIS MANUAL BEFORE OPERATING WITH THE EQUIPMENT.

Please bear in mind the following:

This manual is inseparable from the Nahita digital leukocyte counter, so it should be available for all the users of this equipment.

You should carefully handle the equipment avoiding sudden movements, knocks, free fall of heavy / sharp objects on it. Avoid spilling liquids inside the equipment.

Never dismantle the different pieces of the equipment to repair it yourself, since it could produce a defective use of the whole equipment and a loss of the product warranty, as well as injuries on people that handle the counter.

To prevent fire or electric discharges avoid dry or dusty environments. In case it may happen unplug the equipment immediately.

If you have any doubt about setting up, installation or functioning do not hesitate in contacting your wholesaler. You can also tell us any doubts or suggestions you have by contacting Nahita Technical Assistance Department by email to asistencia@auxilab.es or by telephone: +34 807 117 040 (0.30 Euros/min).

This equipment is protected under the Warranties and consumer goods regulation (10/2003).

Overhaul is not covered by the equipment warranty.

Operations made by non-qualified staff will automatically produce a loss of the warranty.

Neither fuses nor accessories (including their loss), are covered by the product's warranty. The warranty neither covers piece's deterioration due to the course of time.

Please make sure you keep the invoice, either for having the right to claim or asking for warranty coverage. In case you have to send the equipment to Nahita Technical Assistance Department you should enclose the original invoice or a copy as guarantee.

Manufacturer reserves the right to modify or improve the manual or equipment.

ATTENTION!! IF EQUIPMENTS ARE NOT PROPERLY CLEAN AND DISINFECTED THEY WOULD NOT BE ALLOWED TO REPAIR BY OUR TECHNICAL SERVICE.